


UNI*FUSE™ INFUSION CATHETERS

Now available
4F UNI*FUSE™



*Patented slit design
makes all the
difference...*


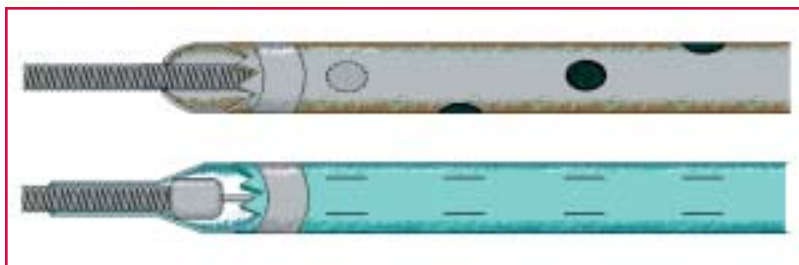


ANGIODYNAMICS®

INCORPORATED

The UNI*FUSE™ advantage... *Everything you like in the 5F is now available in the 4F size!*

Occluding wire prevents kinking

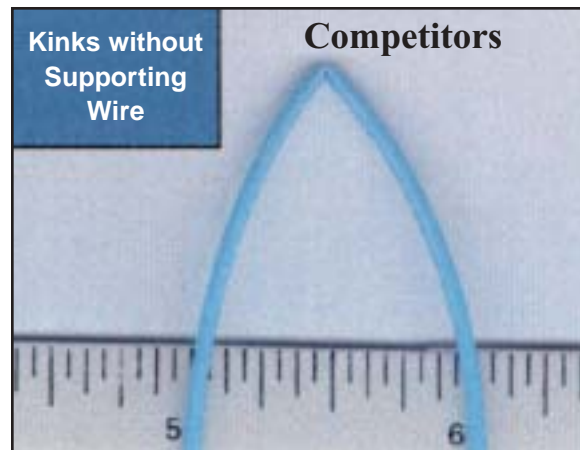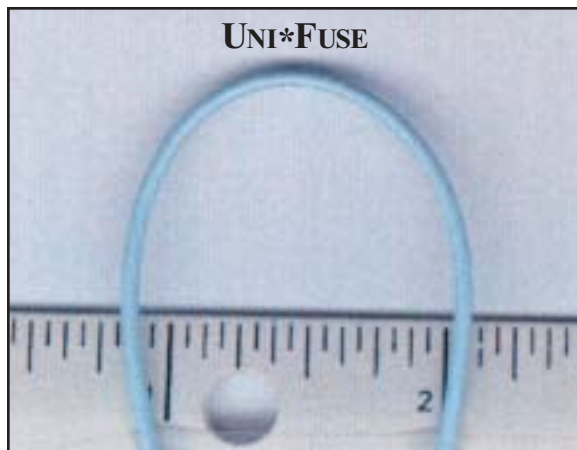

- Supports catheter over the bifurcation



- 0.035" guidewire: gives you more pushability and mechanical advantages
- Radiopaque markers: for ease of placements

INFUSION CATHETERS


- Flow-Thru Hub™ with self-adjusting occluding wire
- Occluding ball wire occludes end hole allowing pressure responsive slits to open



- Slit design—even flow of distribution
- Allows end hole or side slit fluid distribution

Even distribution...in or out of clot

- Tests have shown the UNI*FUSE Infusion Catheter provides even distribution... even under variable infusion conditions.

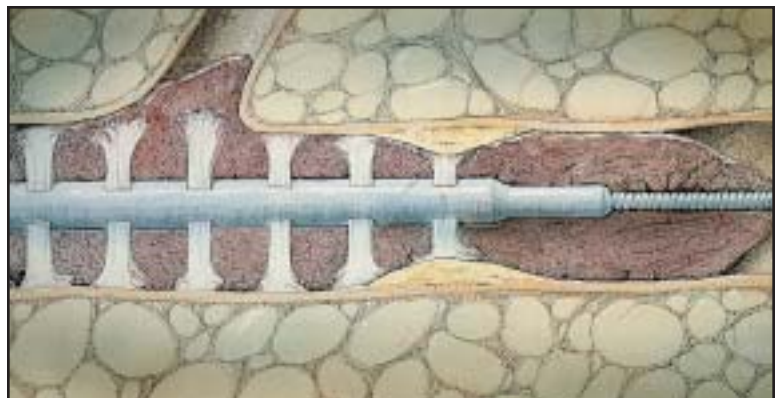


ORDERING INFORMATION

UNI*FUSE™

The only infusion catheter you will ever need...

Flow-Thru Hub™ with self-adjusting occluding wire

- Eliminates Y-adaptor.
- Supports catheter over bifurcation.
- Allows end hole or side slit fluid distribution.
- 4F available in 90cm and 135cm lengths with infusion slit patterns of 5, 10 and 20cm.
- 5F available in 45cm, 90cm and 135cm lengths with infusion slit patterns of 5, 10, 20, 30, 40 and 50cm.

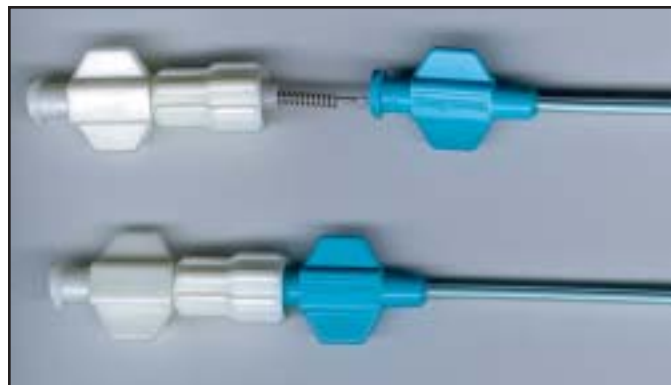

UNI*FUSE™ Infusion Catheter Ordering Information

| Catalog No. | Description | Box |
|-------------|--------------------------------------|-----|
| 12401817 | 4F x 90 cm x 5 cm Infusion Pattern | 1 |
| 12401818 | 4F x 90 cm x 10 cm Infusion Pattern | 1 |
| 12401819 | 4F x 90 cm x 20 cm Infusion Pattern | 1 |
| 12401820 | 4F x 135 cm x 5 cm Infusion Pattern | 1 |
| 12401821 | 4F x 135 cm x 10 cm Infusion Pattern | 1 |
| 12401822 | 4F x 135 cm x 20 cm Infusion Pattern | 1 |
| 12401801 | 5F x 45 cm x 5 cm Infusion Pattern | 1 |
| 12401802 | 5F x 45 cm x 10 cm Infusion Pattern | 1 |
| 12401803 | 5F x 45 cm x 15 cm Infusion Pattern | 1 |
| 12401804 | 5F x 45 cm x 20 cm Infusion Pattern | 1 |
| 12401805 | 5F x 90 cm x 5 cm Infusion Pattern | 1 |
| 12401806 | 5F x 90 cm x 10 cm Infusion Pattern | 1 |
| 12401807 | 5F x 90 cm x 20 cm Infusion Pattern | 1 |
| 12401808 | 5F x 90 cm x 30 cm Infusion Pattern | 1 |
| 12401809 | 5F x 90 cm x 40 cm Infusion Pattern | 1 |
| 12401810 | 5F x 90 cm x 50 cm Infusion Pattern | 1 |
| 12401811 | 5F x 135 cm x 5 cm Infusion Pattern | 1 |
| 12401812 | 5F x 135 cm x 10 cm Infusion Pattern | 1 |
| 12401813 | 5F x 135 cm x 20 cm Infusion Pattern | 1 |
| 12401814 | 5F x 135 cm x 30 cm Infusion Pattern | 1 |
| 12401815 | 5F x 135 cm x 40 cm Infusion Pattern | 1 |
| 12401816 | 5F x 135 cm x 50 cm Infusion Pattern | 1 |

MEDICAL SPECIALTIES
AUSTRALIA PTY LIMITED
54 GIBBES ST, CHATSWOOD NSW 2067
PO BOX 764, WILLOUGHBY, NSW 2068
PHONE 02 9417 7955 FAX 02 9417 5779
E-Mail: custserv@msa.com.au
www.msa.com.au



ANGIODYNAMICS®

INCORPORATED

ANGIODYNAMICS, INC. 603 Queensbury Ave., Queensbury, NY 12804, U.S.A. Phone: 518-798-1215
Uni*Fuse and Flow-Thru Hub are trademarks of AngioDynamics, Inc. U.S. Patents: 6,283,950; 5,250,034.

Fax: 518-798-1360

Toll-Free USA: 1-800-772-6446

www.angiodynamics.com