


THD PressProbe

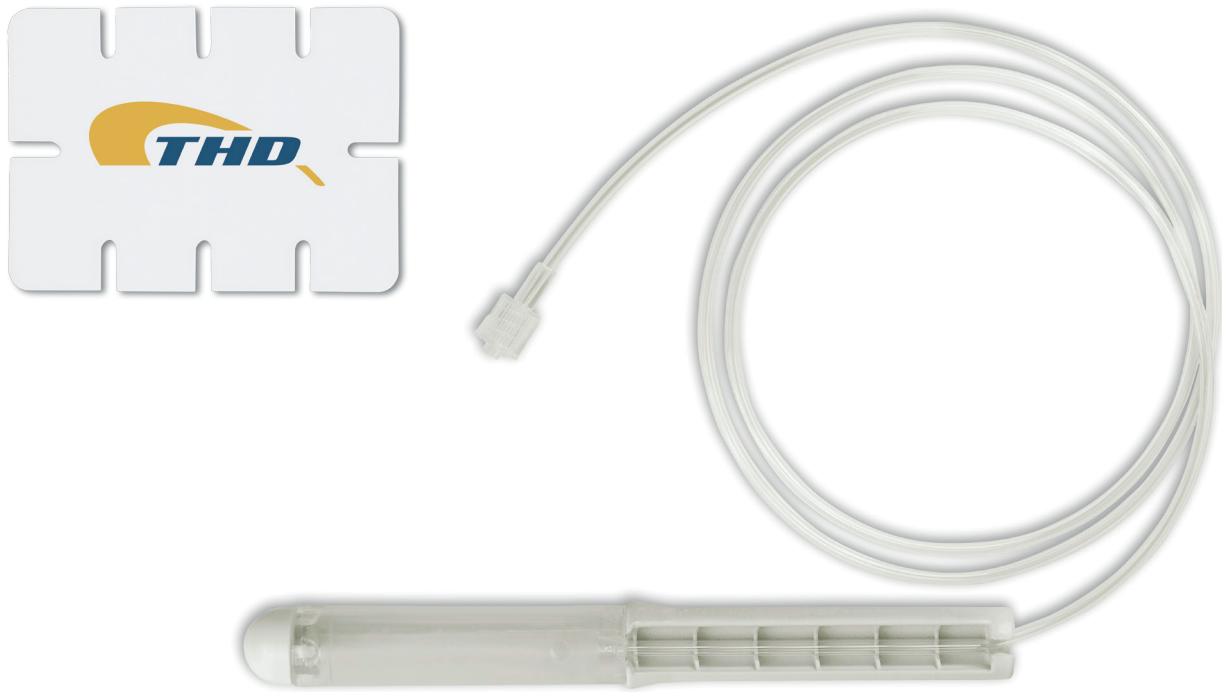

THD PressProbe is the innovative disposable active probe for the execution of the THD anal manometry exam.

The probe allows for the measurement of the real medium sphincter tone pressure in the different functional activations (basal, squeeze and strain tones).

THD PressProbe is meant to be used exclusively in combination with THD Anopress.

There is no need of calibrating the probe, because of the self-calibrating routine before the start of the procedure.

The probe has an ergonomic grip and a completely atraumatic profile, which ensures the maximum comfort for both doctor and patient during the examination and it has been designed to reproduce a normal digital examination.

Each pack contains a proprietary Probe Card, to enable the recording of the examination.

THD PressProbe

Probe Characteristics	
Probe Diameter	16 mm
Probe length	165 mm (including the handle)
Length of measuring space	50 mm
Length of tubing	1000 mm
Probe Weight	23 g
Connection to Anopress	Luer lock

Table of Materials	
Handle + Body	ABS
Obturator + Membrane	PU
Tubing	PU
Luer Lock Connector	MABS

THD PressProbe is available as:

Ref.	Description	Specifications
800150	THD PressProbe (BOX 10)*	Pack of 10 PressProbes in a single box
800150	THD PressProbe (ENV 50)	Pack of 50 PressProbes in a single plastic envelope

*Available for trial only