






A Low-Profile Design. A High Level of Patient Comfort.

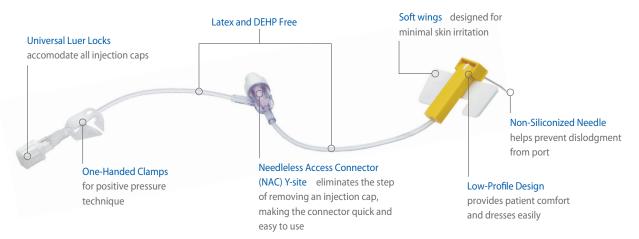

Safe Has Never Been This Easy









Raise needle cover ninety degrees.

Description 19 Ga. x .75" 19 Ga. x 1" 20 Ga. x .50" 20 Ga. x .75" 20 Ga. x 1" 20 Ga. x 1.5" 22 Ga. x .50" 22 Ga. x .75" 22 Ga. x 1"

Hold base firmly.

□ H787LG221005

Lift wings straight up.

>> Select models Power-injectable with the Smart Port (refer to product list for codes)

| UPN | Power Injectable |
|----------------|------------------|
| □ H787LG19755 | Yes |
| □ H787LG191005 | Yes |
| □ H787LG20505 | Yes |
| □ H787LG20755 | Yes |
| □ H787LG201005 | Yes |
| □ H787LG201505 | Yes |
| □ H787LG22505 | No |
| □ H787LG22755 | No |
| | |

| Description | UPN | Power Injectable |
|-----------------------------|-------------------|------------------|
| 19 Ga. x .75" w/ NAC Y-site | □ H787LG1975NY5 | No |
| 19 Ga. x 1" w/ NAC Y-site | □ H787LC19100NX5 | No |
| 20 Ga. x .50" w/ NAC Y-site | □ H787LG2050 (1)3 | Мо |
| 20 Ga. x .75" w/ NAC Y-site | □ H781 LS2075NY5 | No |
| 20 Ga. x 1" w/ NAC Y-site | DH787LG20100NY5 | No |
| 20 Ga. x 1.5" w/ NAC Y-site | □ H787LG20150NY5 | No |
| 22 Ga. x .50" w/ NAC / site | ☐ H787LG2250NY5 | No |
| 22 Ga. x .75" w/ NAC Y-site | ☐ H787LG2275NY5 | No |
| 22 Ga. x 1" w/ NAC 🗸 site | ☐ H787LG22100NY5 | No |

PHYSICIAN'S SIGNATURE

* CT

IMPORTANT RISK INFORMATION

INDICATION FOR USE: The AngioDynamics LifeGuard Safety Infusion Set is designed to access implanted vascular ports to administer fluids and/or to withdraw blood. The LifeGuard Safety Infusion Set is a 90° non-coring needle composed of wings, tubing, a universal Luer connector, optional Needleless Y-site, tubing clamp(s), needlestick guard and needle trap that facilitates the safe removal of the needle by encapsulating the needle during vascular port de-accessing (to help prevent needlestick injuries).

CONTRAINDICATIONS: Do not use with patients showing signs or symptoms of infection at the insertion site or who have exhibited prior intolerance to the materials of

WARNINGS AND PRECAUTIONS: Review all fluids to be used with this product for any incompatibility with plastics and plastic adhesives. A 10 ml syringe or larger is recommended. Clamp infusion set when not in use. Device is intended for single use only. It is non-pyrogenic

and non-toxic. Inner contents are sterile unless package is opened or damaged. Do not resterilize. Carefully read and follow all instructions prior to use of the device. Failure to do so may result in patient complications.

CAUTION: Federal (USA) law restricts this device to sale by or on the order of a physician.



USA > 14 Plaza Drive, Latham, NY 12110 > tel: 800-772-6446 or 518-798-1215 > fax: 518-798-1360 International > Haaksbergweg 75 (Margriettoren), 1101 BR, Amsterdam Z-O > The Netherlands tel: +31 (0)20 753 2949 > fax: +31 (0)20 753 2939

www.angiodynamics.com





SPECIALTIES AUSTRALASIA

MEDICAL

54 Gibbes St, Chatswood 2067 NSW, Australia PO Box 764, Willoughby 2068 NSW, Australia

2 02 9417 7955

www.msa.com.au

■ 02 9417 5779

© 2014 AngioDynamics, Inc. MLC 391 US Rev C